

Pulsox-300i Pulse Oximeter

*“The performance, comfort & memory
you demand from a pulse oximeter”*



*Fits comfortably
on the wrist*




► Pulsox-300i Pulse Oximeter Overview

The **Pulsox-300i** offers the performance, comfort and memory you demand from a pulse oximeter. This unit is excellent for overnight monitoring or anywhere the need for data memorization is required.

PART NUMBERS

Pulsox-300i.....	R124P07
<i>(without probe)</i>	
Pulsox-300i.....	R204P17-002
<i>(with Soft Tip probe)</i>	
Pulsox-300i with finger probe.....	R204P17
<i>(with finger probe)</i>	
Pulsox-300i.....	R204P17-001
<i>(with spot check probe)</i>	

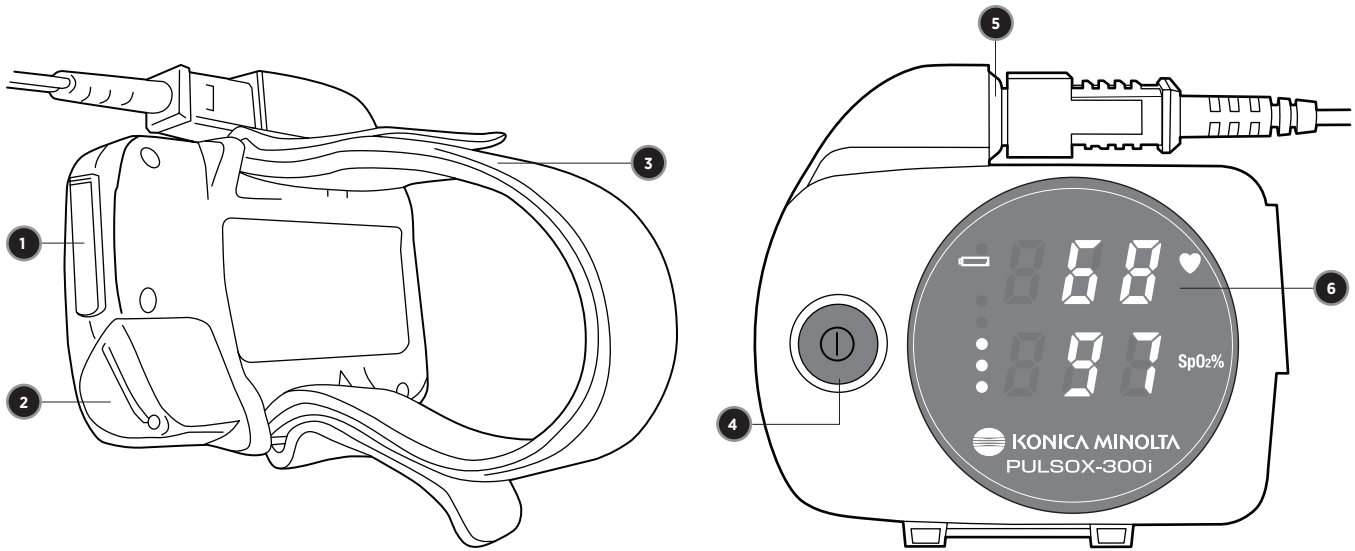

1. External Output Terminal

For exporting data to a computer

2. Battery Cover

Up to 30 hrs with 1 AAA-size alkaline battery

3. Wristband

Fits comfortably on the wrist

4. On/Off Button

Hold down for 3 seconds

5. Probe Connector

Connects to a wide range of SpO₂ probes

6. Large, Back Lit LCD

Provides clarity in low light environments

Specifications

Measuring method:	Dual-wavelength pulse-type oximeter
Measuring range:	SpO ₂ : 0 to 100%, Pulse rate: 30 to 230 bpm
Accuracy:	SpO ₂ : ± 2% (70 to 100%, 1 S.D.), Pulse rate: ± 2 bpm (30 to 100 range) or ± 2% of value (100 to 230 range)
Warning functions:	SpO ₂ warning (flashing display), Battery indication, Error message indication
Memory function:	Battery indication, Error message indication (Up to 300 hours at one-second intervals)
Temperature/humidity range:	0 to 40°C (32 to 104°F); 30 to 85% relative humidity with no condensation
Power:	1 AAA-size battery (approx. 30 hrs with alkaline battery, continuous use)
Dimensions (W x H x D):	68 x 58 x 15 mm (2.68" x 2.28" x 0.59")
Weight:	56 g (2 oz.) (including battery/excluding wrist band and probes)
Standard accessories:	Wristband, 1 AAA-size alkaline battery
Optional accessories:	Soft Tip Probe, Finger Clip Probe, Spot Check Probe, Extension cable, USB adapter, PROFOX software