
Venti Plus Test Lungs

SENSING

ANALYSIS

DELIVERY

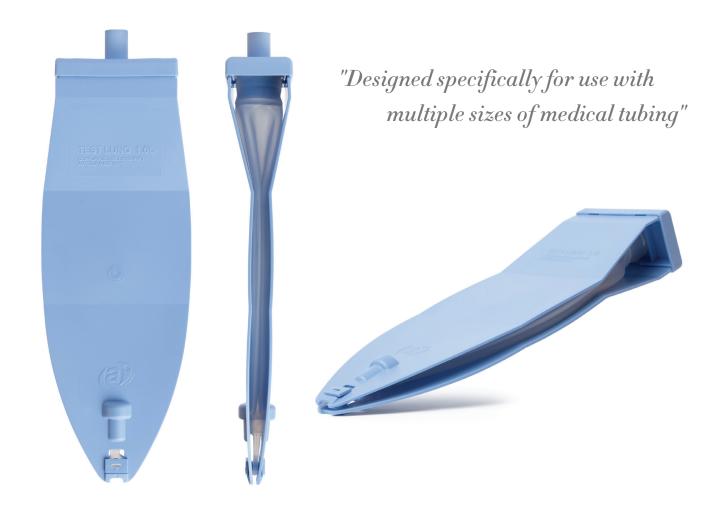

VENTI PLUS TEST LUNG OVERVIEW

The **Venti Plus Test Lung** was designed for fast and easy ventilator circuit testing prior to clinical use. It comes complete with a latex-free silicone bag that simulates either an adult lung (1 liter models), or an infant lung (0.5 liter models). Both models include a turn and lock design that allows for easy removal of the bag for cleaning and decontamination. Every lung is calibrated for accuracy and validated for resistance and compliance within the application range.

PART NUMBERS

Specifications

VENTI PLUS 1.0L (ADULT)

Compliance 0.02 L/cmH₂O

5 cmH₂O/L/Second Resistance

Volume 0-600 ml (with 1L bag)

Weight 230 grams

Dimensions 300 x 100 x 43mm

Warranty 12 months

VENTI PLUS 0.5L (PEDIATRIC)

0.01 L/cmH₂O Compliance

5 cmH₂O/L/Second Resistance

Volume 0-400 ml (with 0.5L bag)

Weight 140 grams

240 x 100 x 43mm **Dimensions**

12 months Warranty

