


# Universal Equipment Stand

SENSING

ANALYSIS

DELIVERY

*"Accommodates most  
CPAP & BiPAP units"*

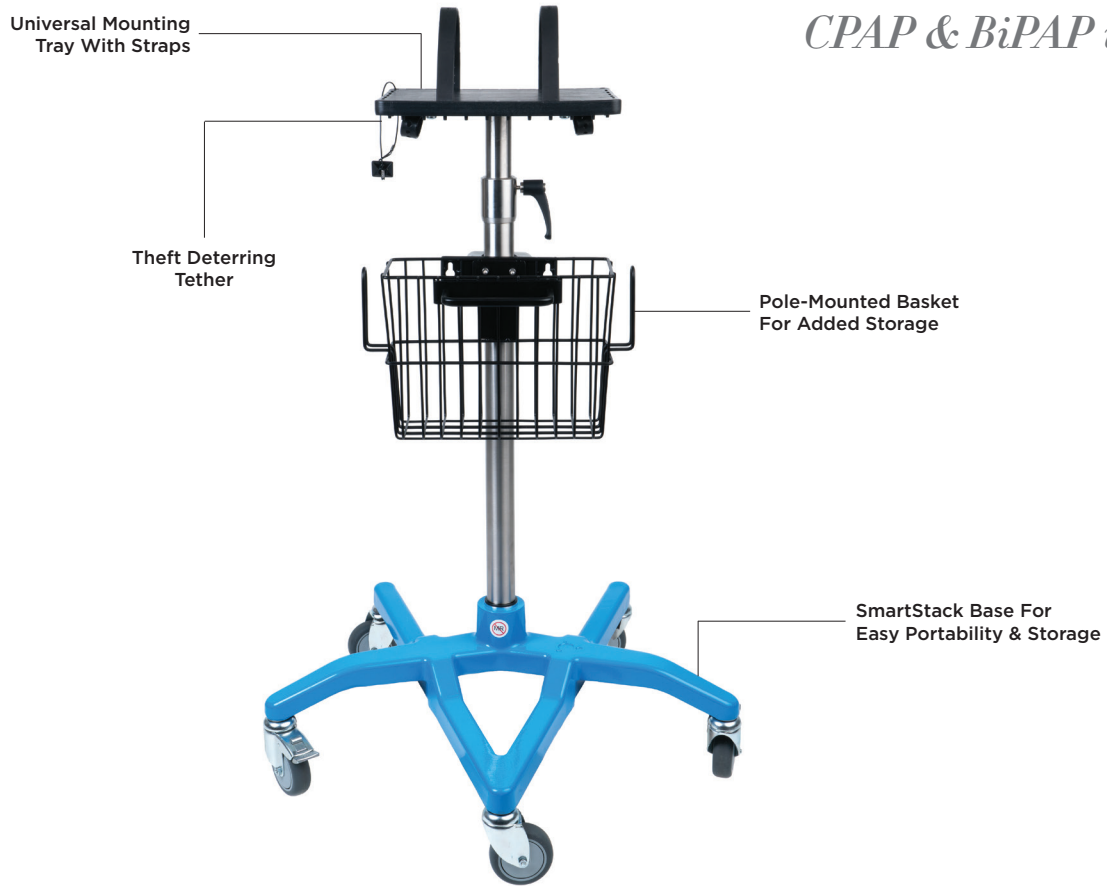

## A UNIVERSAL EQUIPMENT STAND


The **SmartStack** universal equipment stand comes equipped with a heavy-duty, universal mounting plate, which is designed to accommodate most CPAP and BiPAP units. Its specially-designed cast aluminum base provides extra support and durability. Its patented interlocking design nests neatly with other **SmartStack** bases for ease of portability and storage. The **SmartStack** equipment stand also includes a pole-mounted basket for added storage.

## PART NUMBERS

Universal Equipment Stand . . . . . R100P86

\*Patent# 8,196,874 B2

\*Japan# 5214740e

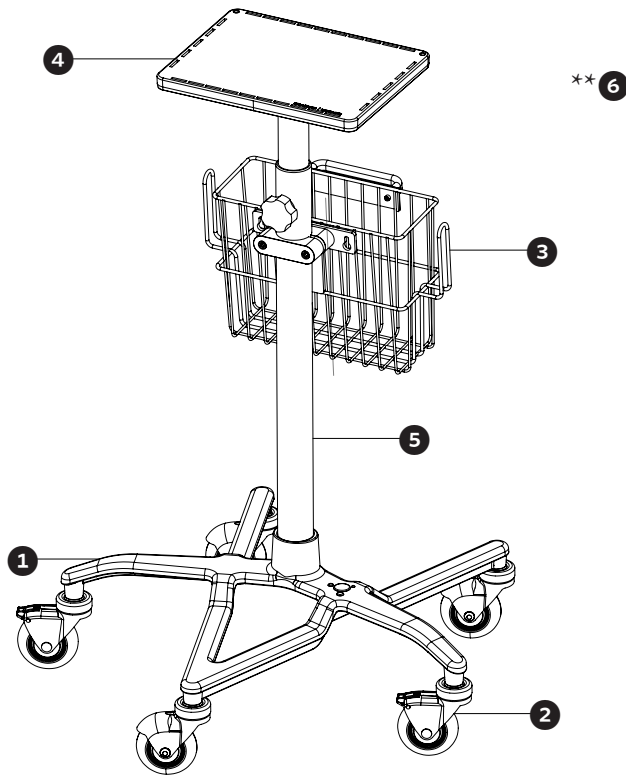

**1. SmartStack Base**

*Durable Aluminum Base*

**2. Locking Casters**

*Large, Smooth-Rolling Casters With 2 Wheel Locks*

**3. Wire Basket**

*Guidance Handle & Cord Hooks*

**4. Heavy Duty Universal Mounting Tray (8.5" x 11")**

*Securely Mount CPAP & BiPAP Units to Stand With Straps (Straps Included)*

**5. Chrome Plated Pole**

*1.75" Stainless Steel Pole*

**6. \*\*Mounting Hardware**

*Straps & Other Misc. Mounting Methods Included (listed below at bottom of page)*

**Specifications**

|                          |   |
|--------------------------|---|
| Collapsed Height:        | 34" (Measured from floor to mounting plate)                       |
| Extended Height:         | 45" (Measured from the floor to the mounting plate)               |
| Casters:                 | Neoprene Tread, Chrome Steel, 3" QTY: 5, 2 Locking, 3 Non-Locking |
| Weight Capacity          | 30 lbs total  |
| Wheel Base Diameter      | 22"   |
| Nesting Base             | 8" center to center   |
| Weight                   | 20 lbs  |
| Shipping Weight          | 12 lbs (Pole, basket, plate) 12 lbs (Base and casters)            |
| Pole Diameter            | 1.75" (Bottom portion) 1.5" (Upper Portion)                       |
| Base Material            | Powder Coated Aluminum  |
| Pole Material            | Stainless Steel   |
| Restrictions             | Not MRI Compatible  |
| Warranty                 | 48 Months   |
| Mounting Tray Dimensions | 8.5" X 11"  |

*\*\*Universal Equipment Stand comes equipped with Universal Mounting Tray, screws to mount tray to SmartStack Equipment Stand and means to secure equipment to tray: elastomeric straps, velcro, lanyard and adhesive pad, and a drill template on which mounting locations of the equipment may be marked and transferred to the tray for drilling holes for equipment mounting screws.*